



瀏覽更多
禮品即入
www.giftu.com.hk



Inspire your gift idea
www.giftu.com.hk

MATRIX PROMOTION LIMITED

Ray Leung
852 – 90495142
mark@matrixgroup.com.hk
www.giftu.com.hk

Modern and stylish reusable mask

XDDESIGN



Price 網購

Home | Price 網購 | 免費註冊 | 瀏覽紀錄 | 項目比較 | 項目分類 | 帳號 | 二手 | Hot | 網購 | 購物中心 | 語言 | App下載

搜索 網購商品

熱度: 預購電視 | Nintendo Switch | Switch | BOSE | 按摩椅 | 洗衣機

網購首頁

- 本週人氣熱賣
- 精選蔬菜
- Mobo+開食
- 限時精選
- 夏季外出必備
- 最新上架
- \$200有找
- 影音產品
- 數碼產品
- 電腦產品
- 家庭電器
- 家品傢俬

RHYTHM 雙葉 手提風扇 \$299 \$198

XD Design 納米纖維 口罩套裝 \$299 \$188

全港免費 送貨上門

本週熱賣

- 日本 Vdada 智能防蟻電飯煲
- XD Design 納米纖維口罩套裝
- Iris Ohyama 超強全方位靜音
- Sony WF-1000XM3 真無線降噪
- 日本 COLUMIDA alex 菌治恐龍




XD Design Masks are designed for extra comfort



- 3-dimensional chin structure, a nose bridge that endures for thousands of bends
- Adjustable straps - Free the pain on the ears
- Slim & Lightweight

1. OUTSIDE FABRIC - polyester knitted fabric treated with ViralOff®
2. 1ST FILTER LAYER - PP non-woven
3. 2ND FILTER LAYER - Hyproof® Nanofiber microporous membrane
4. 3RD FILTER LAYER - PP non-woven
5. INNER FABRIC - polyester knitted fabric

XD Design Protective Masks - Treated with ViralOff®

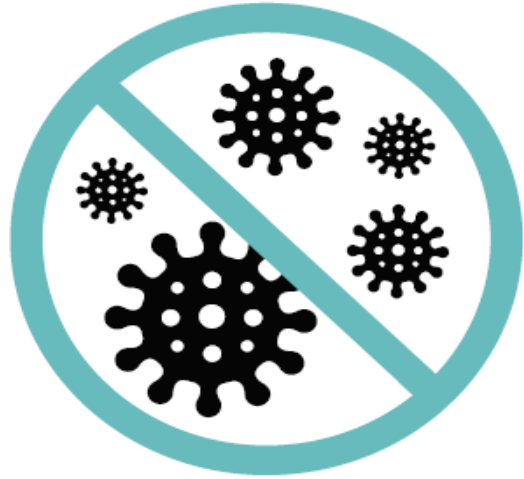

99% REDUCTION OF
VIRUSES IN 2 HOURS

ISO18184 tests against viruses (like Influenza, Corona, Hepatitis)



ViralOff
POLYGIENE TECHNOLOGY



ViralOff
POLYGIENE TECHNOLOGY

ViralOff® is a technology that stops viral activity through interaction with key proteins. It shows a 99% reduction of many different viruses, on a treated material, within a 2 hours timespan. Shown by ISO18184:2019, ATCC VR-1679 (H3N2) standard antiviral test.

The active ingredient present in ViralOff® is a reaction mass of titanium dioxide and silver chloride.

Due to ViralOff treatment the XD Design mask can be used during all day. It's not necessary to wash the mask on every use.

(However It is not recommended to wear the mask longer than 2-4 hours for each use).

The mask (with filter removed) is machine washable on 60 degrees and it can be washed depending on the frequency of use .

The ViralOff® is a durable and always-on treatment. For best performance and sustainability, wash less and only when needed.



XDDESIGN x Viraloff Protection Reusable Mask

Key features:

Kill germs and virus – 99% (SGS testing report)

Hyproof® nano fibers filters x 5 included , **used up to 1000 hrs**

Carrying pouch included Washable

Custom logo available



Product Size :15*11cm
Material: Nanofibers
Imprint: Pad Printing
Packing: pouch





XD Mask - Testing Report - P265.87 breathability (mask + filter)

Mask -test

Nano fiber Filter
BFE > 96%, PFE > 96%

Replacable filters



SGS
 Test Report SLS2025269630301TX Date: June 17, 2020 Page 2 of 2

Medical Face Masks-Requirements and Test Methods (EN 14885:2016+AC:2019)

Clause 5.2.2 Breathability (EN 14885:2016+AC:2019 Annex C, Flow rate 8 l/min)

	1#	2#	3#	4#	5#
Differential pressure ΔP (Pa/cm²)	71	70	77	69	72

Remark: Performance Requirement: Type I<40 Pa/cm², Type II<40 Pa/cm², Type III<60 Pa/cm²

Sample Photo



The statement of conformity in this test report is only based on measured values by the laboratory and does not take their uncertainties into consideration.

End of Report

SGS **ILAC-MRA** **CNAS**
 Test Report SLS2025269630301TX Date: June 17, 2020 Page 1 of 2

The following sample(s) was/were submitted and identified on behalf of the client as:
 Sample Description : XD Design protective mask set

SGS Internal Ref No. : SHH-L2006520745MD
 Sample Color : Grey
 Country or Destination : EU
 Item No. : P265.87 XD DESIGN Protective Mask Set

Test Performed : Selected test(s) as requested by applicant

Sample Receiving Date : Jun 03, 2020
 Testing Period : Jun 04, 2020 - Jun 17, 2020
 Test Result(s) :

Signed for and on behalf of
 SGS-CSTC Standards Technical Services (Shanghai) Co., Ltd Testing Center

Sara Guo
 Sara Guo (Account Executive)

CTI CHINA CERTIFICATION & INSPECTION GROUP SHENZHEN CO., LTD.
 Report No. 4478420050710 Date: 06/04/2020

Test Requested :

Test Items	Bacterial Filtration Efficiency (BFE), Filtration efficiency
Testing basis	YY 0469-2011 Surgical mask

Test Results:

Test item(s)	Test method(s)	Limit	Measured Values	Judge
Bacterial Filtration Efficiency (BFE) (%)	YY 0469-2011 Appendix B (Bacterial species: Staphylococcus aureus ATCC 6539; Mean particle size: 3.0 μm; Flow rate: 28.3 L/min; Test area: 40cm²; Positive Control: 1.94 × 10 ⁶ CFU; Negative Control: <1CFU)	≥ 95	95.6 95.2 95.8	Pass
	Filtration efficiency (%)	≥ 30 (Non oily particles)	99.0 99.5 99.2	Pass

Note: Pass = Up to standard
 Fail = Not up to standard
 "-" = Not regulated
 For the entrusted sample test, the technical responsibilities are undertaken for the submitted samples only.
 The test results of this report are only for internal quality control reference.

Approved by: *[Signature]*
 检测技术有限公司

SGS
 Test Report SLS2025269619501TX Date: June 17, 2020 Page 2 of 2

Medical Face Masks-Requirements and Test Methods (EN 14885:2016+AC:2019)

Clause 5.2.2 Bacterial filtration efficiency (BFE) (EN 14885:2016+AC:2019 Annex C)

(BFE), %	1#	2#	3#	4#	5#
	99.9	99.5	99.9	99.9	99.9

Remark: Performance Requirement: Type I:≥95%, Type II:≥96%, Type III:≥96%

Sample Photo



The statement of conformity in this test report is only based on measured values by the laboratory and does not take their uncertainties into consideration.

End of Report

HS CODE : 6307900010.

The carton size is 61.5 x 39 x 32.5 cm, 180pcs /per carton, net weight 5.40kg, gross weight 6.40kg, unite product set is 30g, box sizes is 15.5 x 12.5 x 2 cm.

XDDESIGN x Viraloff Protection Reusable Mask

Product highlights

1. The technology: the mask itself is **made with Viral off ® fabric**, as briefly explained in the presentation the fabric is treated with a material **that kills up to 99% of the bacteria and viruses** that are caught in the fabric. It's therefor self cleaning in a sense.
2. 1 mask per set (**available in 3 different colors**) with **5 filters** and a pouch to carry the mask and filter when not using them. The pouch is made of corrugated carton and is simple in use.
3. The filters are made of **Hyproof® nano fibers that help with proper filtration** but more importantly are washable up to 10x. Each filter can be used up to 200 hrs, so with the package you' ll have a **1000 hrs of usage**. The product is in no way of medical grade but adheres to Afnor Spec S76-001



Replacable filters



Portable pouch



Pre-order price:

- **500pcs – HK\$62.3/pcs**
- **1000pcs - HK\$58.7/pcs**
- **2000pcs – HK\$55.9/pcs**



Kwun Tong Station



APM

HOW MING STREET

HOI YUEN ROAD

SHING YIP

HUNG TO ROAD

HING YIP STREET

Welcome visit our showroom

“MAKE YOUR OWN GIFT !”

MATRIX PROMOTION LIMITED

(852) 3620 3717 joe@matrixgroup.com.hk
10/F, The Rays, 71 Hung To Street, Kwun Tong, Kowloon, H.K.

**GiftU
Showroom**

